

使用方法說明

安立復美達®

持續性藥效肌肉注射用懸浮劑


預充填
注射筒

400 mg

每個月注射一次




套組內容： 安立復美達注射筒

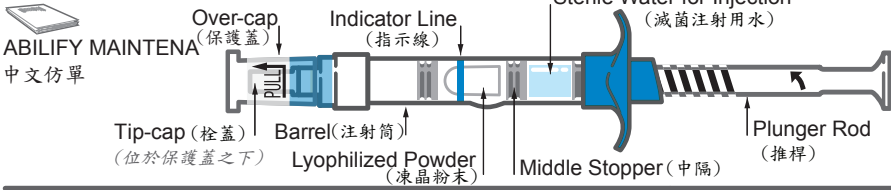

- ▲ 1吋針頭
(供非肥胖，三角肌使用)
- ▲ 1.5吋針頭
(供非肥胖，臀肌使用
或肥胖，三角肌使用)
- ▲ 2吋針頭
(供肥胖，臀肌使用)

1 釋出稀釋液(diluent)


推動及旋轉
推桿(plunger rod)
直到推桿完全停止。

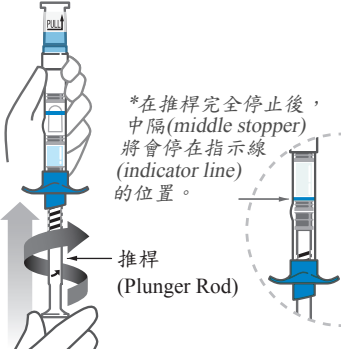


*在推桿完全停止後，
中隔(middle stopper)
將會停在指示線
(indicator line)
的位置。

推桿
(Plunger Rod)

2 垂直地用力地振搖注射筒20秒，直到藥液呈現
均勻的乳白色。

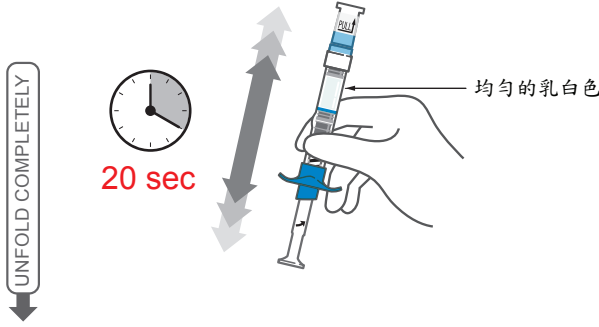


3 扭轉(Twist)及拉掉
注射針筒保護蓋
(Over-cap)及
栓蓋(Tip-cap)

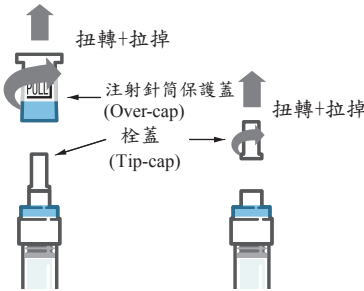


4 依照體型及注射部位，選擇適當的注射針。

體型	注射部位	針頭大小
非肥胖	三角肌(Deltoid)	1 inch (23G)
	臀肌(Gluteus)	1.5 inch (22G)
肥胖	三角肌(Deltoid)	1.5 inch (22G)
	臀肌(Gluteus)	2 inch (21G)

5 抓住針頭套(needle cap)，以順時鐘方向轉動，
確保注射針牢牢地固定在注射針筒上。然後以
垂直向上的方式將針頭套自針頭上拉掉。

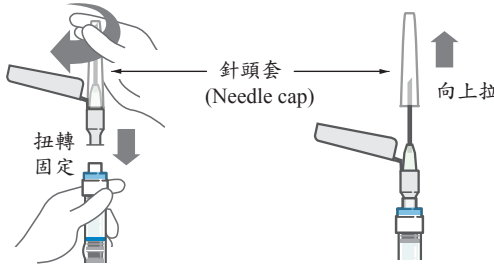


6 垂直地握住注射針筒，然後慢慢地推動推桿
(plunger rod)以排出空氣。

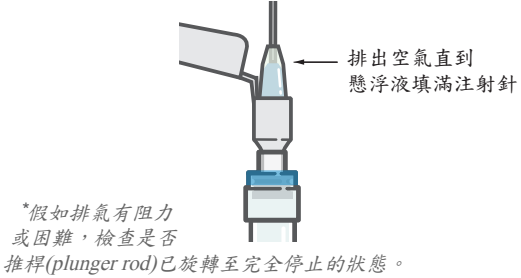


*假如排氣有阻力
或困難，檢查是否
推桿(plunger rod)已旋轉至完全停止的狀態。

7 慢慢地注射ABILIFY MAINTENA至適當的
注射部位。
藥液在再配製後立即注射。

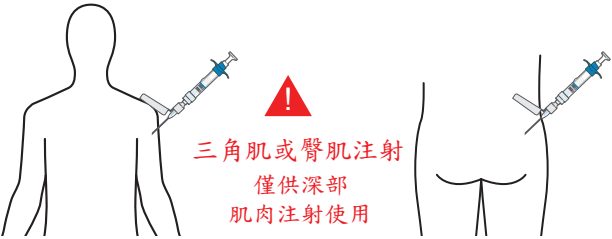

8 確定注射針已套入安全裝置，然後安全地丟棄所有的
套組元件。

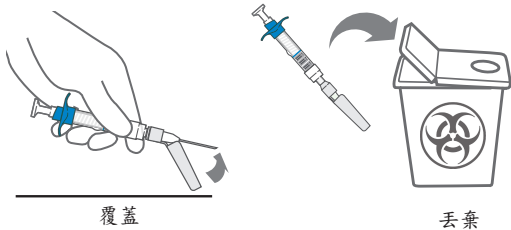


Quick Reference Guide

Abilify Maintena®
(aripiprazole) for extended release injectable suspension


Pre-Filled
Dual Chamber
Syringe

400 mg

Administer once monthly

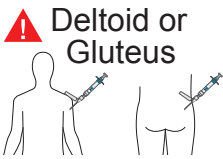


Kit Contents:

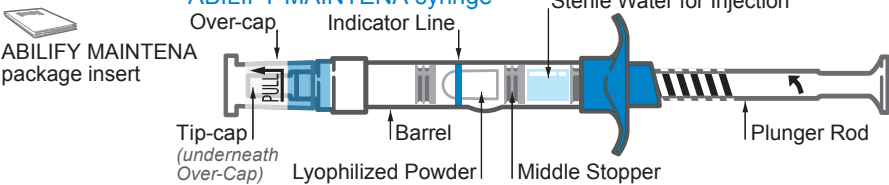

- ▲ 1 inch needle (for non-obese, deltoid)
- ▲ 1.5 inch needle (for non-obese, gluteus or obese, deltoid)
- ▲ 2 inch needle (for obese, gluteus)


1

To release diluent **PUSH and ROTATE** plunger rod **until plunger rod stops** completely.

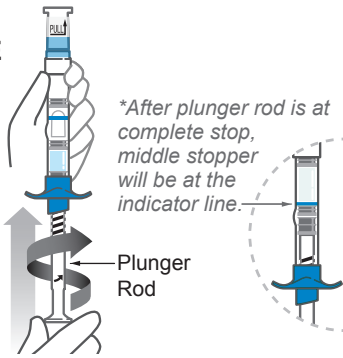


2

Vertically shake syringe vigorously for **20 seconds** until drug is **uniformly milky white**.

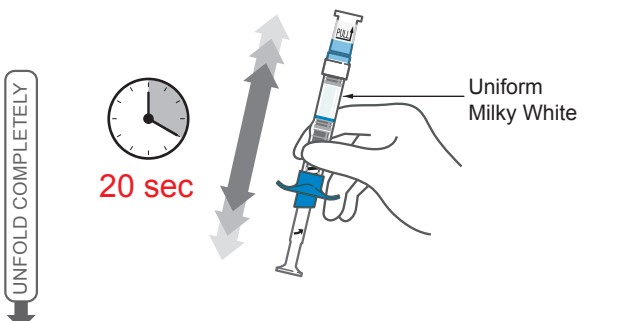

3

Twist and pull off **Over-cap** and **Tip-cap**.

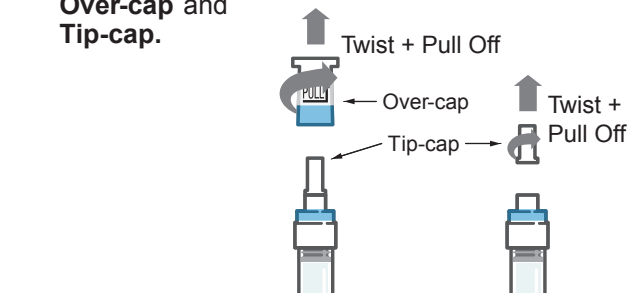


4

▲ Select appropriate needle depending on body type and injection site.

Body Type	Injection Site	Needle Size
Non-obese	Deltoid	1 inch (23G)
	Gluteus	1.5 inch (22G)
Obese	Deltoid	1.5 inch (22G)
	Gluteus	2 inch (21G)


5

Holding needle cap, ensure the needle is firmly seated on the safety device with a push. Twist clockwise until **SNUGLY** fitted. Then **PULL** needle-cap straight up.

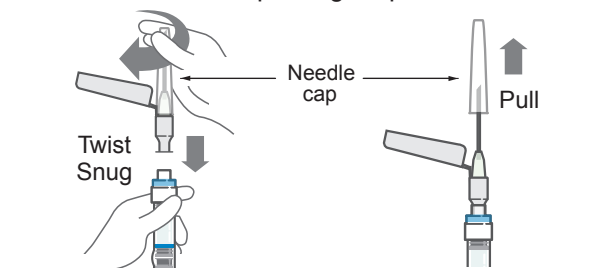


6

Hold syringe upright. **ADVANCE PLUNGER ROD SLOWLY TO EXPEL THE AIR**.

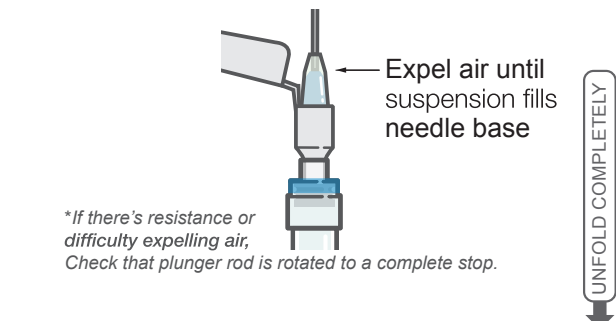

7

Slowly inject ABILIFY MAINTENA into appropriate injection site. Inject immediately after reconstitution.

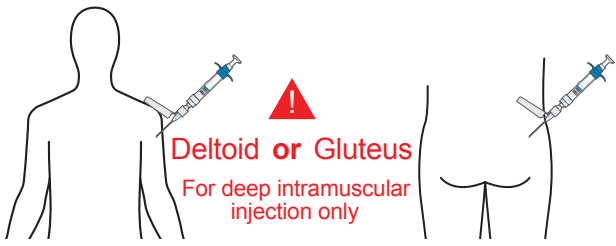

8

Cover the needle and safely discard all kit components.

